

✦ **NEW HANDLE DESIGN:**

Nordent[®]

ONPOINT[®]

RESIN SCALERS AND CURETTES

LIGHTWEIGHT AND COMFORTABLE



New Resin Handle Design

Lightweight handle and ControlRing™ Grip provide excellent comfort and control.

Premium Quality

Nordent tips are made by master craftsmen.

Long-Lasting Sharpness

UltraSteel™ technology results in tips with long-lasting sharpness and strength.

Easy Identification & Customizable

Nordent OnPoint[®] offers **14 silicone ColorRings[®]** that attach to the instrument handle for quick identification.

OnPoint[®]


Explore our wide range of patterns and exclusive Nordent instruments



Nordent UltraSteel™ Technology

UltraSteel technology takes advantage of advancements in metallurgy and cryogenic hardening to produce a micro-structure not previously possible. Scalers made from stainless steel often have large, isolated carbide molecules that can compromise hardness and result in inconsistent wear resistance.

UltraSteel technology features extremely fine carbide particles that are evenly distributed throughout the micro-structure, resulting in tips with long-lasting sharpness and strength.

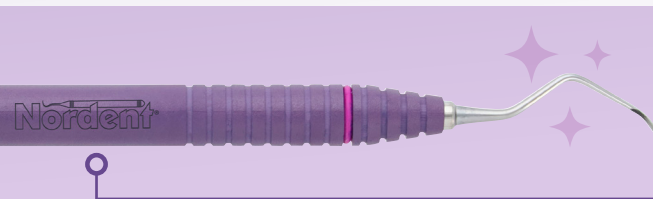

Under **500x Magnification**, these SEM (Scanning Electron Microscope) photos show the **uniform carbide particle distribution** throughout the **micro-structure** of our Nordent UltraSteel technology.


UltraSteel™ Technology
Scalers & Curettes



Leading
Sharpen-Less Scaler



Other Leading Brands
(440A)

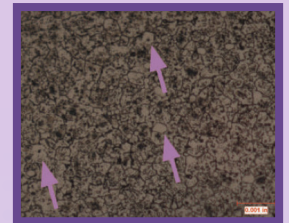

Extensive Selection of Scalers and Curettes



Gracey Curettes

Instrument Tips	Handle Color	Color Ring™
Gracey #1-2	Pink	Gray
Gracey #3-4	Pink	Gray
Gracey #5-6	Pink	Gray
Gracey #7-8	Pink	Gray
Gracey #9-10	Pink	Gray
Gracey #11-12	Pink	Black
Gracey #11-14	Pink	Black
Gracey #12-13	Pink	Black
Gracey #13-14	Pink	Yellow
Gracey #15-16	Pink	Black
Gracey #17-18	Pink	Yellow
Gracey #1-2 Long Reach	Pink	Gray
Gracey #3-4 Long Reach	Pink	Gray
Gracey #5-6 Long Reach	Pink	Gray
Gracey #7-8 Long Reach	Pink	Gray
Gracey #11-12 Long Reach	Pink	Black
Gracey #13-14 Long Reach	Pink	Yellow
Gracey #15-16 Long Reach	Pink	Black

Universal Curettes

Instrument Tips	Handle Color	Color Ring™
Barnhart #1-2	Purple	Purple
Barnhart #5-6	Purple	Purple
Columbia #13-14	Purple	Purple
Columbia #4R-4L	Purple	Purple
McCall #13-14	Purple	Purple
McCall #13S-14S	Purple	Purple
McCall #17-18	Purple	Purple
Univ of CA Rule #3-4	Purple	Purple
Younger-Good #7-8	Purple	Purple
Universal R144	Purple	Purple

Anterior Scalers

Instrument Tips	Handle Color	Color Ring™
Jacquette #1S-Sickle N5 (5-33)	Blue	Blue
Anterior Sickle-Spoon #1	Blue	Blue
Scalette N137 Mini Blade / Long Shank	Blue	Blue
Sickle H5-Hoe #48	Blue	Blue
Sickle H5-H5S Thin	Blue	Blue
Sickle #H5-Spoon #1	Blue	Blue
Anterior Sickle #N67 (H6-H7)	Blue	Blue
Universal Sickle #NMJ (R138)	Blue	Blue
Univ of NEB #128	Blue	Blue

Posterior Scalers

Instrument Tips	Handle Color	Color Ring™
Posterior Sickle #204	Yellow	Yellow
Posterior Sickle #204S	Yellow	Yellow
Posterior Sickle #204SD	Yellow	Yellow
Jacquette #2S-3S (34-35)	Yellow	Yellow
Posterior Sickle #2	Yellow	Yellow
Universal Sickle Offset #4	Yellow	Yellow
Scalette Posterior #N135	Yellow	Yellow
Universal Sickle #N129	Yellow	Yellow
Scaler N130 (MT4-5)	Yellow	Yellow

View our full selection of scaler and curette patterns at NORDENT.COM

© 2026 Nordent and all other marks are either registered or owned trademarks of Nordent Manufacturing, Inc. All rights reserved. Product and specifications are subject to change without notice. DTD #0P7125